

# PowerSpin™ DX Centrifuge


Model: C8724H



- POWERSPIN™ DX CENTRIFUGE
- VARIABLE SPEED 1000-3400 RPM
- 24 PLACE FIXED ANGEL ROTOR
- 60 MIN DIGITAL TIMER
- CAPACITY 24X10 ML TUBES OR 16X15 ML TUBES



The UNICO® PowerSpin™ DX C8700 Series centrifuge is engineered to meet the demanding needs of clinical, industrial, and biological labs. This model comes with an angled rotor designed to hold up to 24 x 10ml tubes. It offers reliable and efficient performance for blood separation and other material processing tasks. This centrifuge features digital variable speed control with a speed range from 500 RPM to 4000 RPM, an electronic timer with a range of 1:00 to 99:59, and a zero RPM safety lid latch that ensures safe operation by preventing the lid from opening while the unit is in motion. The intuitive controls allow users to select from preset programs for blood, urine, and fecal samples, or customize protocols easily using the set, up, down, and enter buttons. The unit's memory function retains the last time and speed settings, ensuring seamless operation for repetitive tasks. With a choice of rotors for standard 3 to 50 ml tubes, and optional adapters for smaller tubes, the C8700 Series is versatile enough for a wide range of laboratory applications. The centrifuge's maintenance-free brushless DC motor delivers quiet, reliable performance, while its rubber suction cup feet ensure stability during operation. Built with safety in mind, the unit will not run unless the lid is securely closed, and it features an automatic time and speed memory reset for convenience. Whether you're working in a clinical or research environment, the UNICO® PowerSpin™ DX C8700 Series offers a reliable and efficient solution for your centrifugation needs.

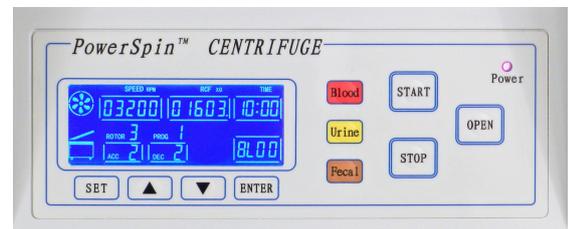

## ADDITIONAL INFORMATION

<b>TOP SELECTION</b>	No
<b>ACCESSORIES</b>	<p><b>Accessories</b></p> <p>C8600-01 6 Place Horizontal Rotor, 2-15ml for DX Centrifuge</p> <p>C8600-02 4 Place Horizontal rotor, 50ml DX Centrifuge</p> <p>C8600-03M 24 Place Fixed Angle Rotor, 2-15ml for DX Centrifuge, Metal (without Tube Shields)</p>
<b>SPECIFICATION</b>	<p><b>Speed (RPM)</b> 500-4000rpm</p> <p><b>Maximum RCF (g)</b> 1848g</p> <p><b>Rotor</b> 24 x 10ml Angled</p> <p><b>Tube Size</b> 2-15ml</p> <p><b>Motor</b> DC</p> <p><b>Power Input</b> 110V</p>
<b>SIZE AND WEIGHT</b>	<p><b>Size</b> 13.75 inch Height (22 inch height with lid open) x 15 inch Deep x 11 inch Width</p> <p><b>Shipping Size</b> 19"X18"X15"</p> <p><b>Weight</b> 31 LBS</p> <p><b>Shipping Weight</b> 36 LBS</p>
<b>WARRANTY</b>	1 years

